

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18


Unison UK-C18

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Donor building block for organic semiconductor ( Halogenated ethylenedioxythiophene analogues )

有機半導体ドナービルディングブロック  
(ハロゲン化エチレンジオキシチオフェン類縁体)

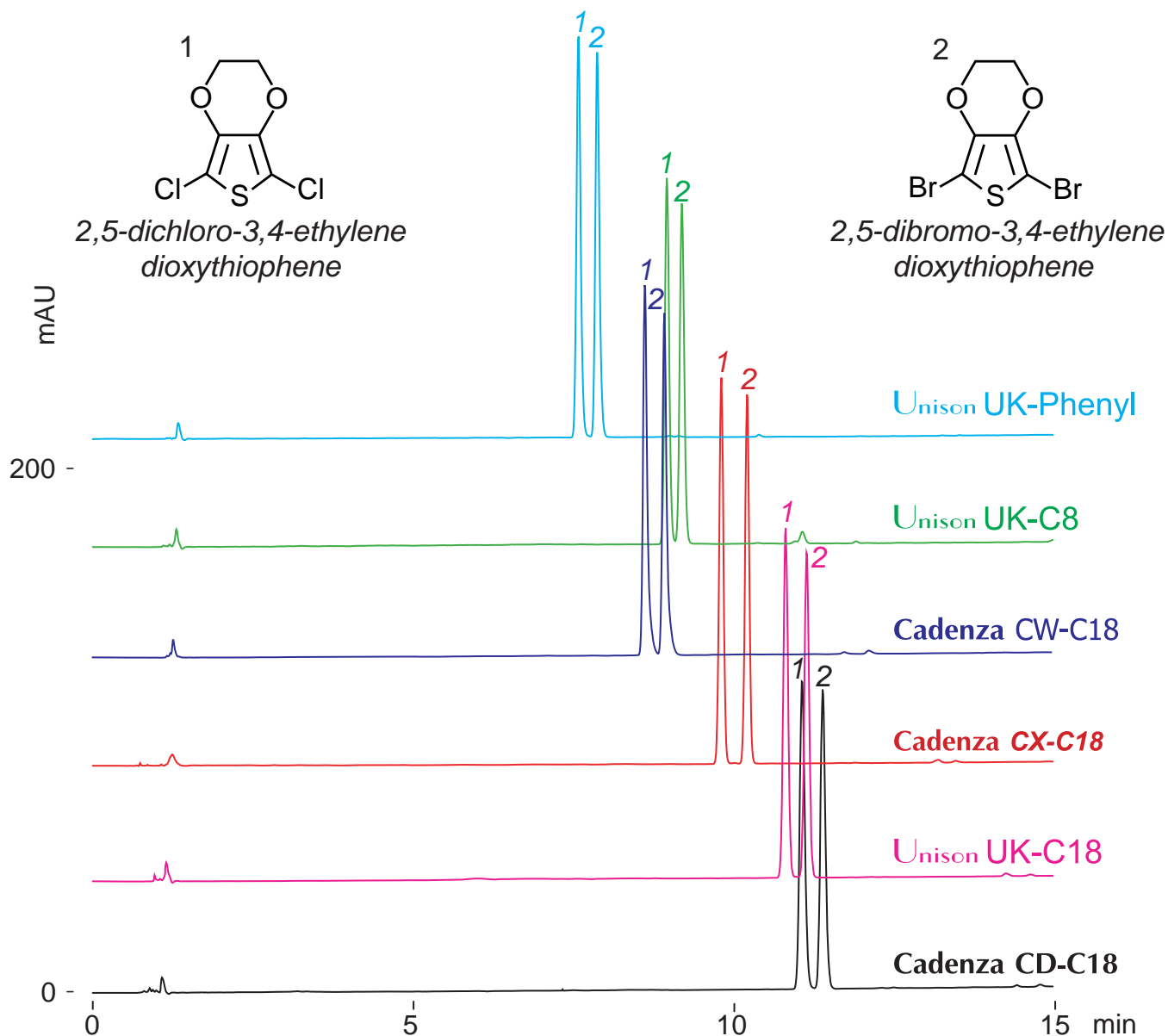

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

20-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.05mg/mL, water / acetonitrile / formic acid = 40 / 60 / 1)