

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18


Unison UK-C18

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Donor building block for organic semiconductor (Dibromocarbazole analogues)

有機半導体ドナービルディングブロック  
(ジブromoカルバゾール類縁体)

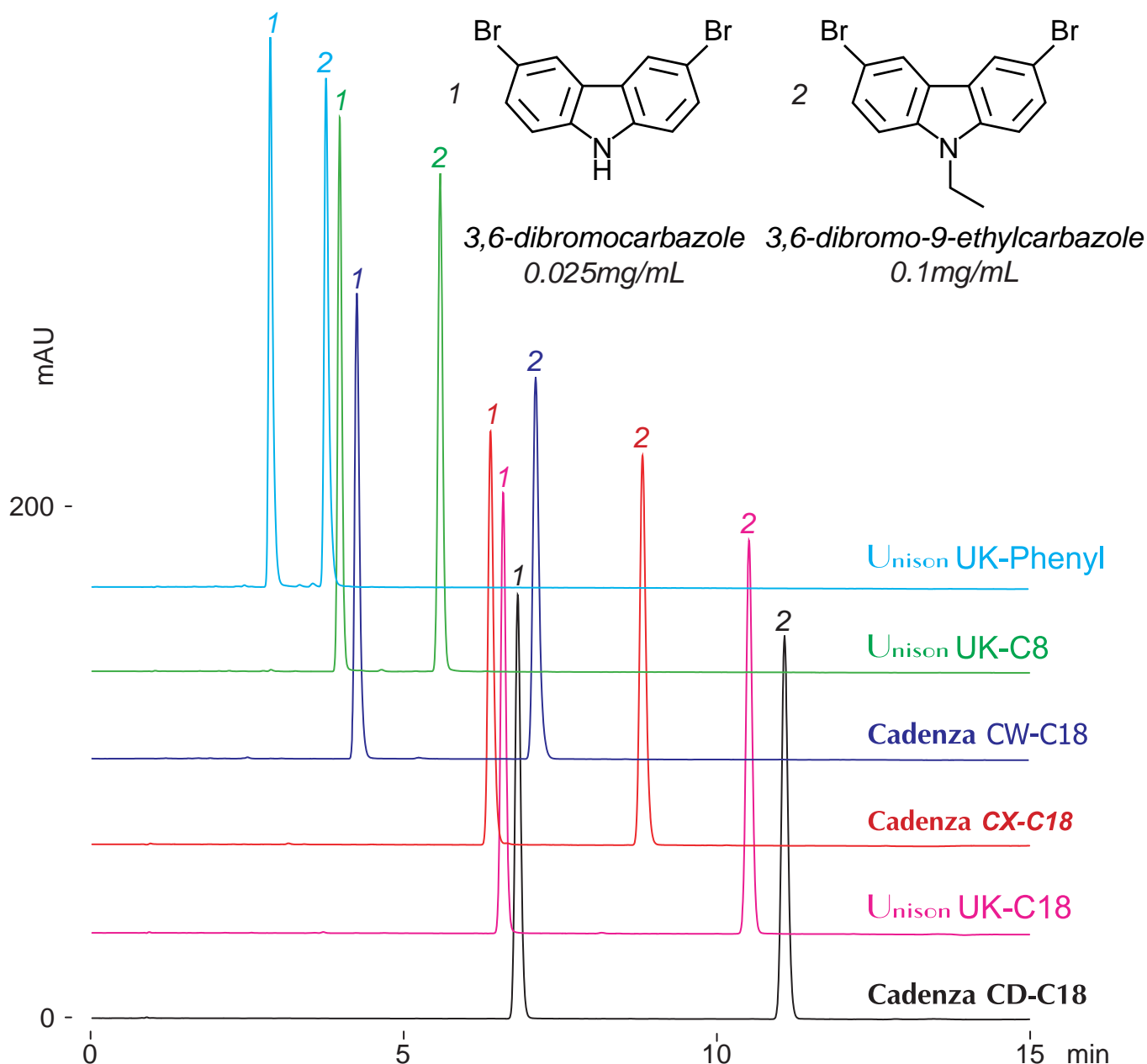

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 50 / 50

B: 100mM ammonium formate / acetonitrile = 10 / 90

30-100%B (0-12min), 100%B (12-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.05-0.1mg/mL, water / acetonitrile / HCOOH = 20 / 80 / 1)