

Unison UK-C18


Cadenza CD-C18

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Gemfibrozil (Hypolipidemic agent)

ゲムフィブロジル (脂質降下薬) の LC-MS 分析

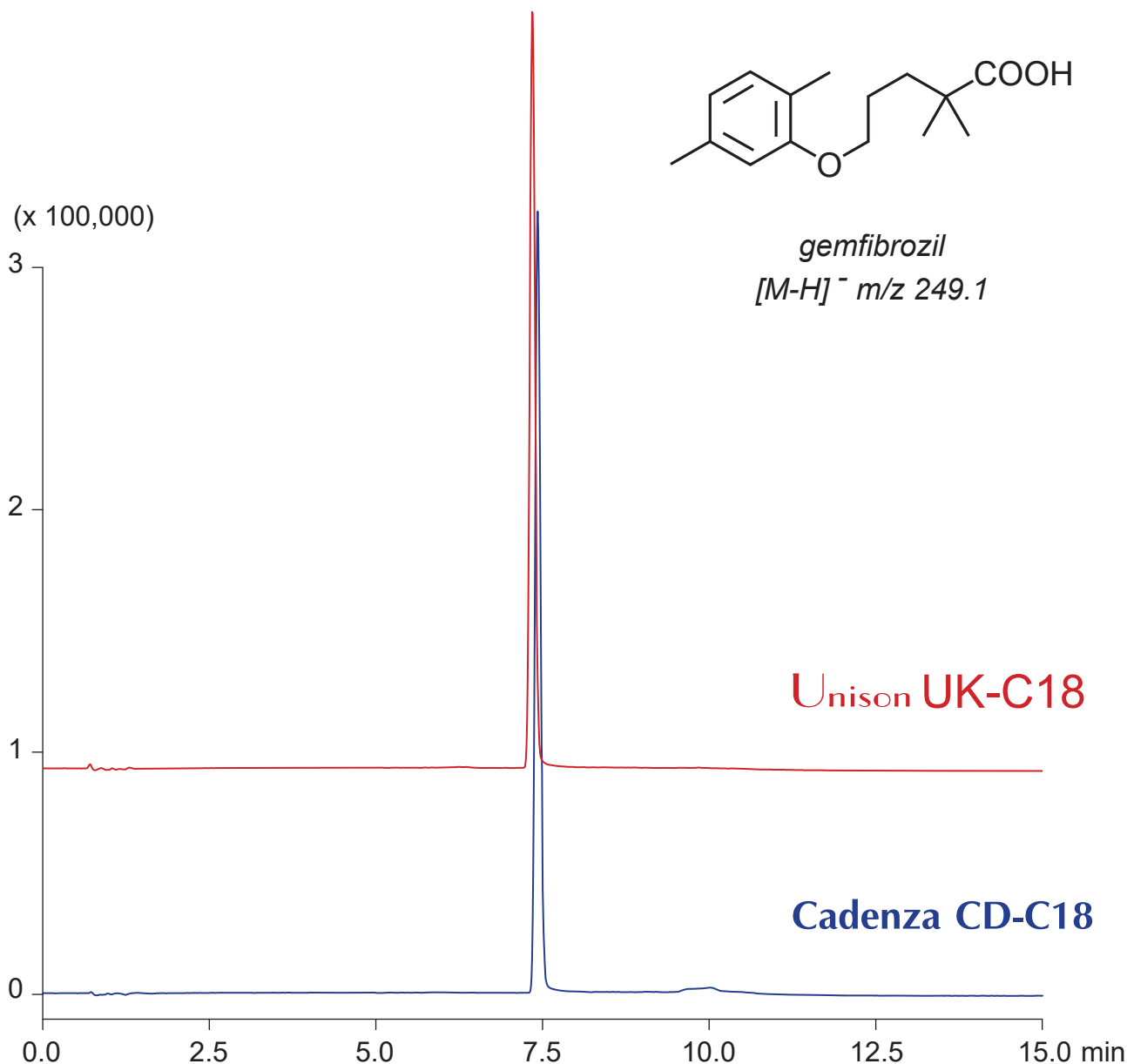

100 x 3 mm

A: water / acetonitrile / acetic acid = 60 / 40 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

0-100%B (0-10min), 100%B (10-15min)

0.4mL/min (7MPa), 37°C, Single Quad. MS (ESI)

5µL (1µg/mL, ethanol / water / acetic acid = 15 / 85 / 1)