


Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION


## Rosuvastatin (Hypolipidemic agent)

ロスバスタチン (高脂血症治療薬)

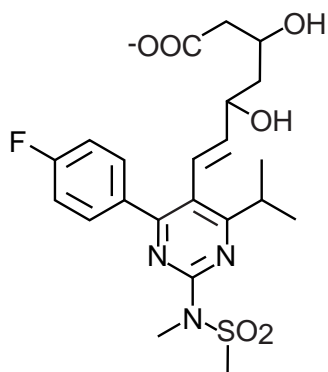

1/2Ca<sup>2+</sup>

rosuvastatin calcium

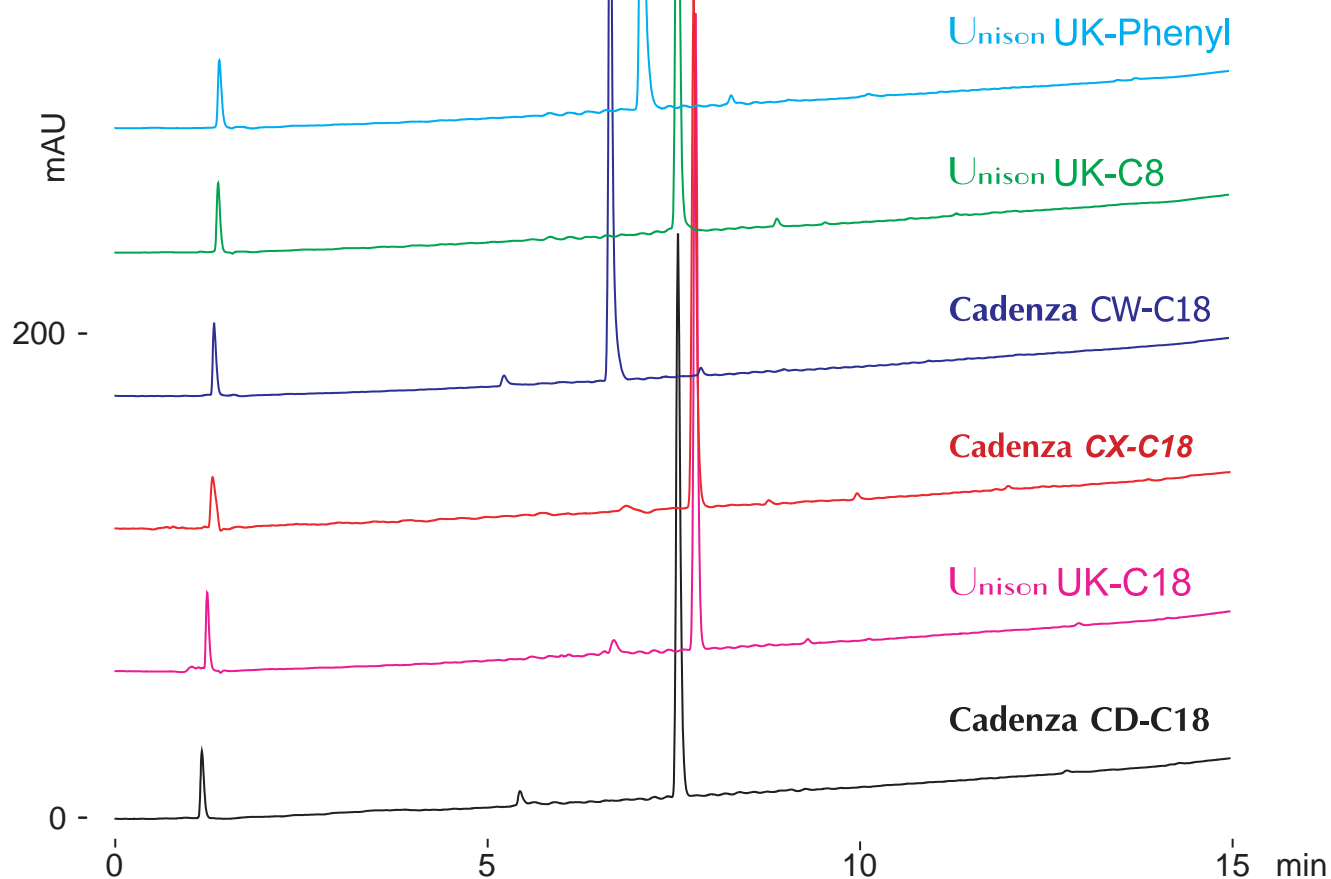

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-90%B (0-15min)

1.0mL/min (10MPa), 37°C, 250nm

10µL (0.1mg/mL, water / acetonitrile / formic acid = 50 / 50 / 1)