


Cadenza CD-C18

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Statins (HMG-CoA reductase inhibitors)

スタチン類 (MG-CoA還元酵素阻害剤) の LC-MS分析

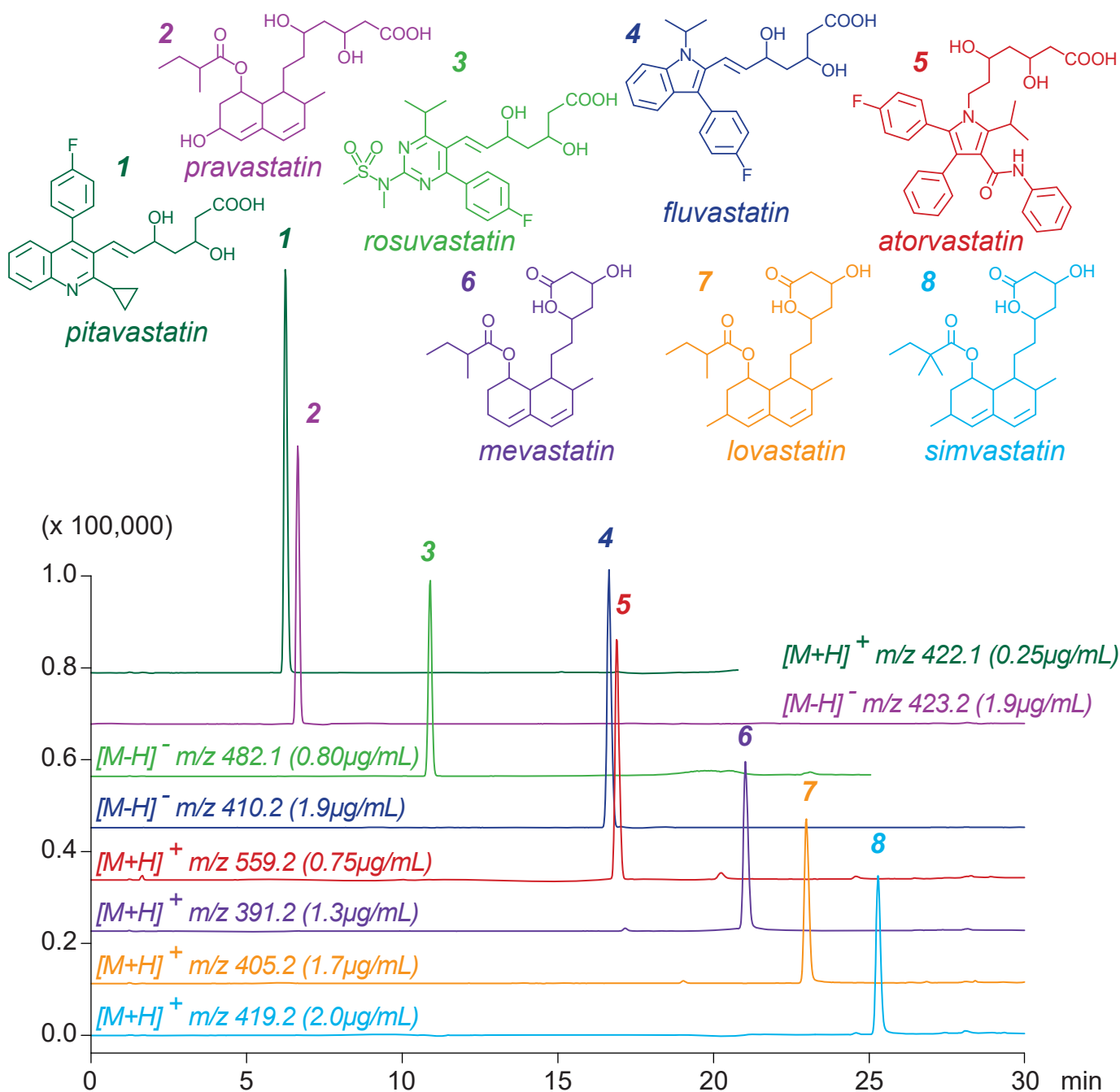

Cadenza CD-C18, 150 x 3 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

25-85%B (0-25min), 100%B (25-30min)

0.4mL/min (10MPa), 37°C, Single Quad. MS (ESI)

 10 μL (0.25-2.0 $\mu\text{g/mL}$, 100mM AcONH₄ / EtOH = 50 / 50)