


Unison UK-Phenyl
 Unison UK-C8
 Unison UK-C18
 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Enzalutamide (Antineoplastic) and synthetic intermediate
 エンザルタミド (抗悪性腫瘍薬) と合成中間体

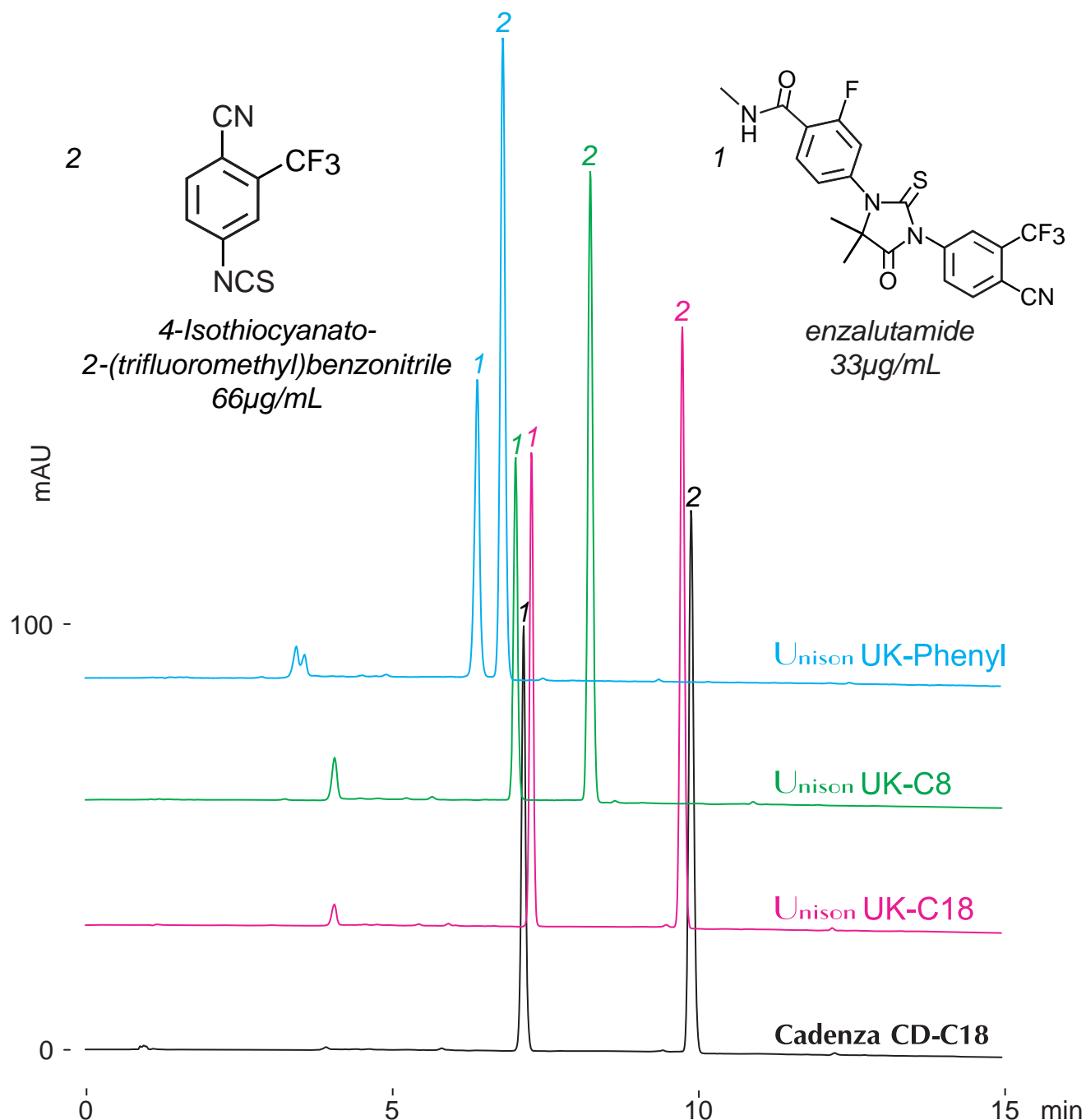

100 x 4.6 mm

A: 100mM ammonium formate = 10 / 90

B: 100mM ammonium formate / acetonitrile = 10 / 90

30-100%B (0-15min)

1.0mL/min (11MPa), 37°C, 260nm

10 μ L (33-66 μ g/mL, acetonitrile)