

Unison UK-C8


Unison UK-C18

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Tris (2,2'-bipyridyl) ruthenium (II) (Dye for dye sensitized solar cell)

トリス (2,2'-ビピリジル) ルテニウム (II) (色素増感太陽電池用色素) の LC-MS 分析

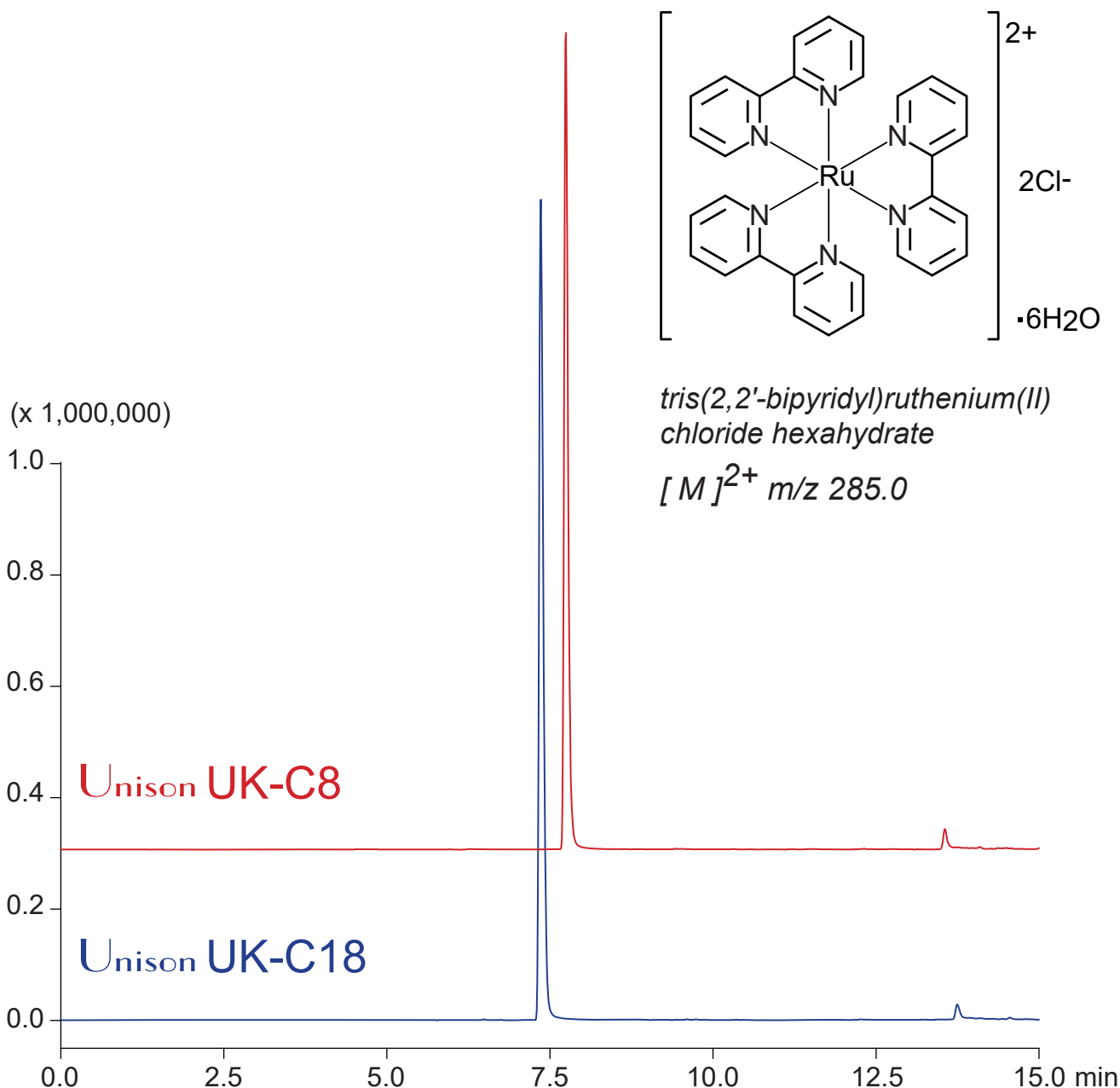

100 x 3 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0%B (0-1min), 0-60%B (1-12min), 100%B (12-15min)

0.4mL/min (7MPa), 37°C, Single Quad. MS (ESI)

10µL (50ng/mL, 100mM ammonium formate)