


Unison UK-C1  
 Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

**Coumarin 153 (Organic light-emitting diode material)**

クマリン153 (有機エレクトロルミネッセンス材料)

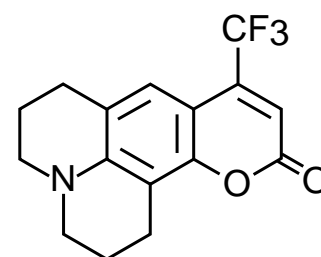

coumarin 153

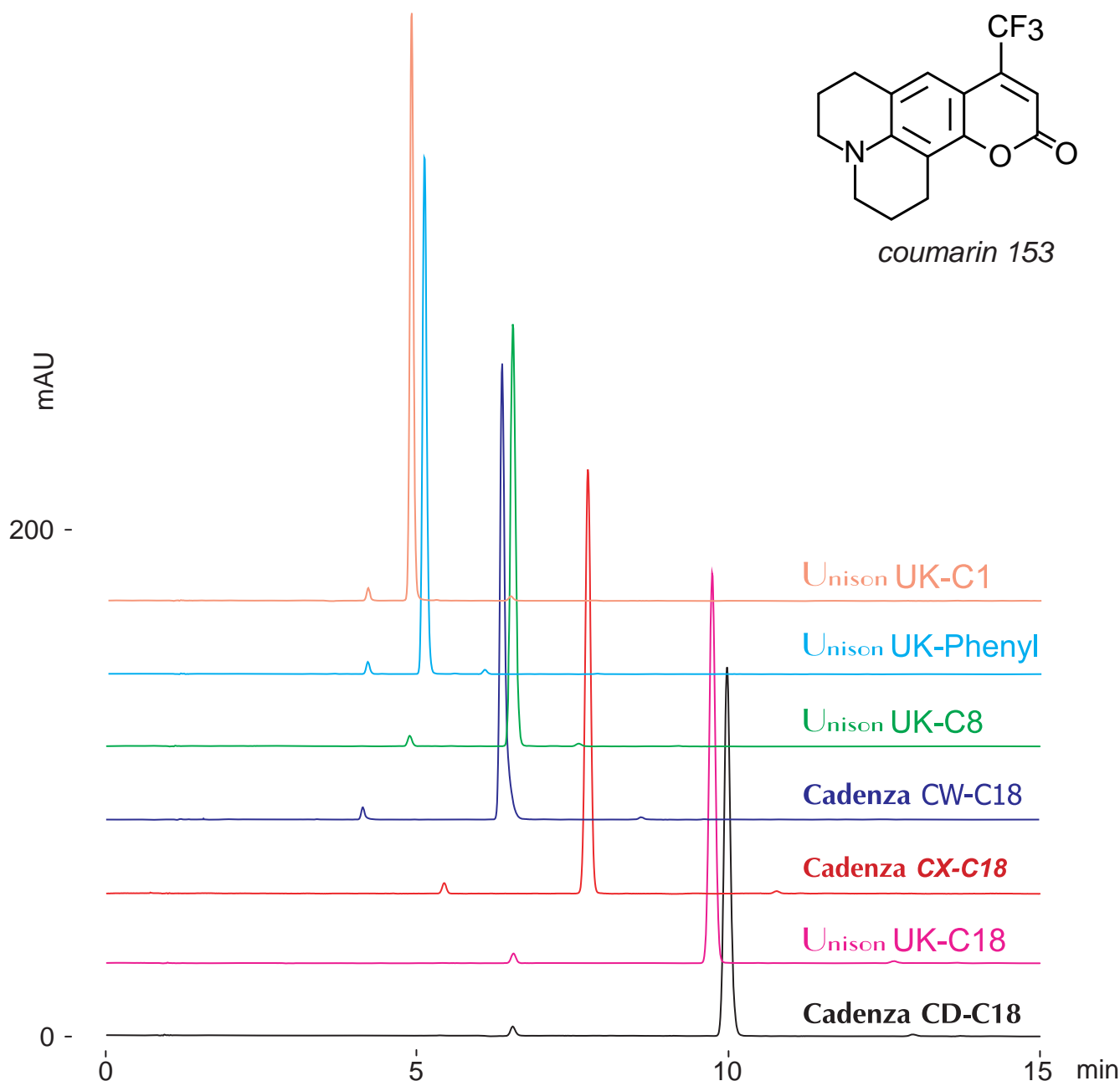

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 410nm

10µL (75µg/mL, water / acetonitrile / formic acid = 35 / 65 / 1)