


Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

1,4-bis [4- (N,N-diphenylamino) styryl] benzene (Organic light-emitting diode material)

1,4-ビス [4- (N,N-ジフェニルアミノ) スチリル] ベンゼン
 (有機エレクトロルミネセンス材料)

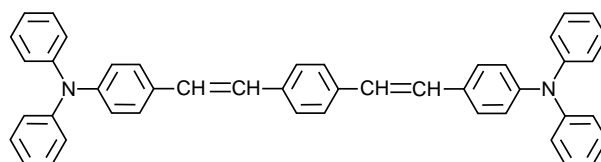

1,4-bis [4- (N,N-diphenylamino) styryl] benzene

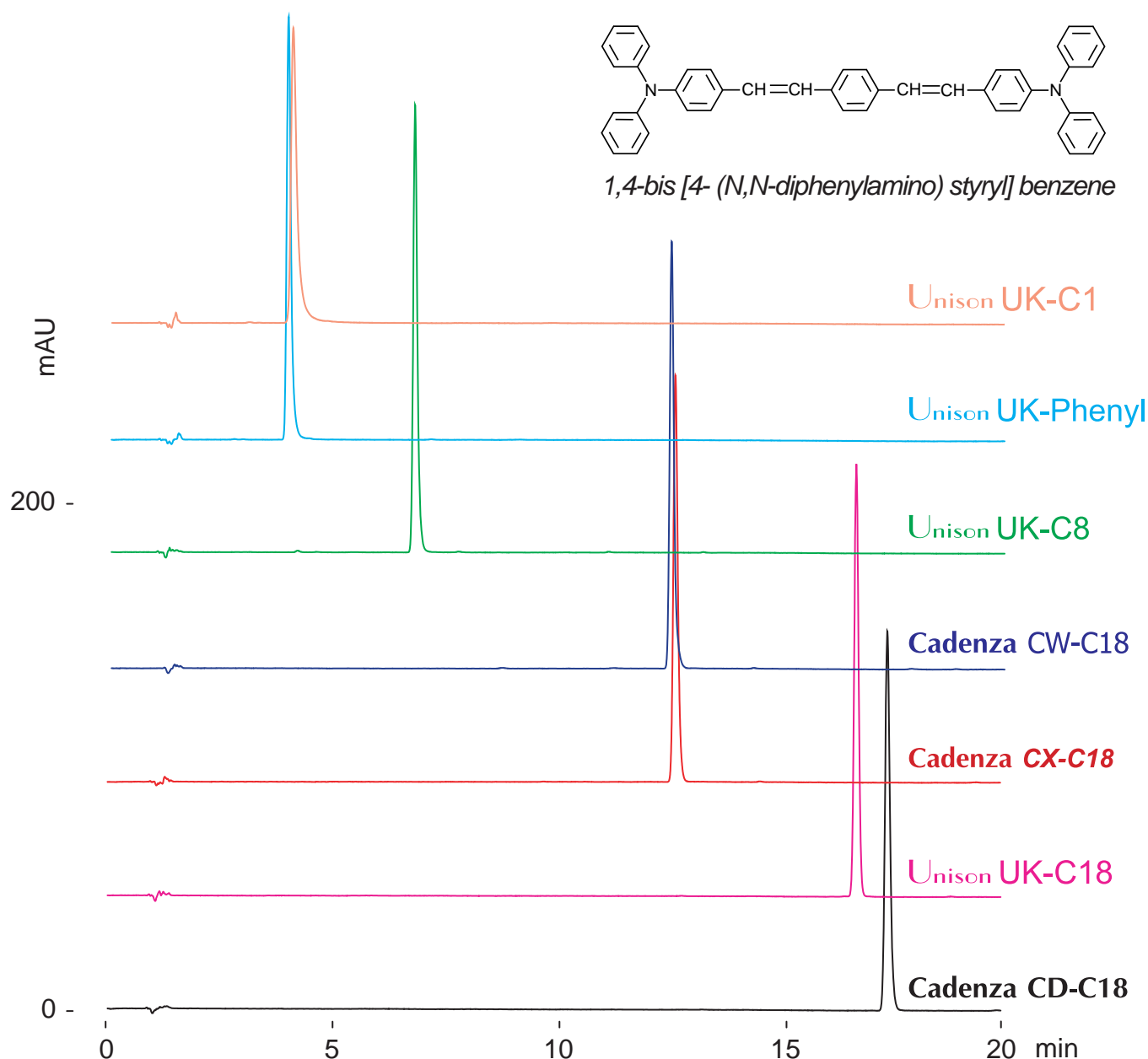

100 x 4.6 mm

A: 100mM HCOONH₄ / acetonitrile = 50 / 50

B: acetonitrile

60-100%B (0-15min), 100%B (15-20min)

1.0mL/min (5MPa), 37°C, 260nm

10μL (0.033mg/mL, DMSO)