


Unison UK-C18  
 Unison UK-C8  
 Unison UK-C1  
 Unison UK-Phenyl

150 x 4.6 mm


HPLC COLUMN APPLICATION

### Sclareol, (+)-sclareolide (Sesquiterpenes)


スクラレオール, (+)-スクラレオリド (セスキテルペン類) の分析

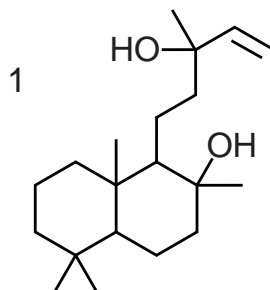

sclareol  
(3mg/mL)

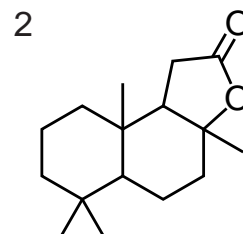

(+)-sclareolide  
(17mg/mL)

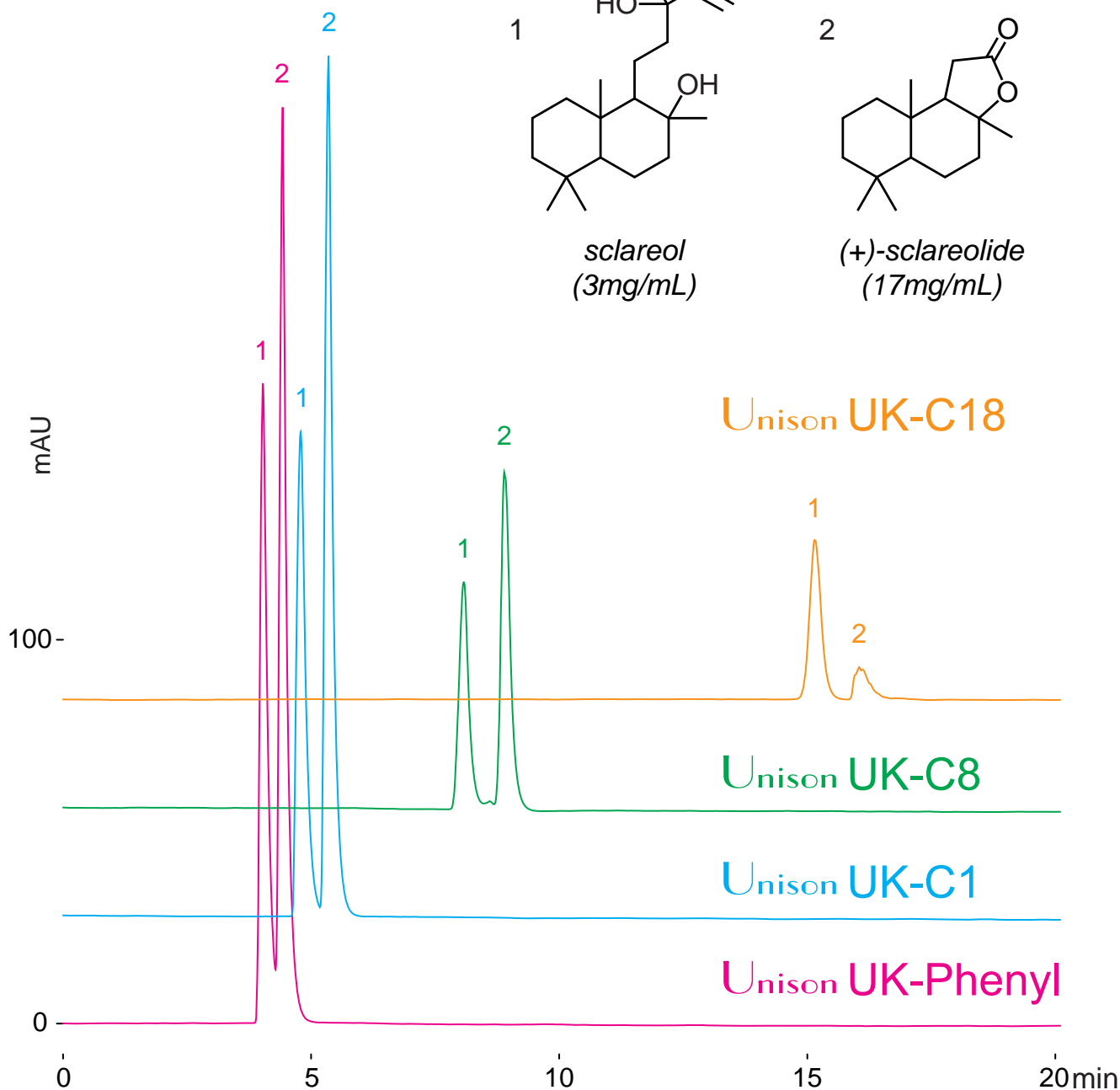

150 x 4.6 mm

acetonitrile / water / acetic acid = 55 / 45 / 0.1

1.0mL/min (10-12MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 70°C)

10µL (3-17mg/mL, tetrahydrofuran)