


Unison UK-C8

250 x 4.6 mm


HPLC COLUMN APPLICATION

1,4-bis [2- (9-ethylcarbazol-3-yl) vinyl] benzene  
(Organic light-emitting diode material)

1,4-ビス [ 2-( 9-エチルカルバゾール -3-イル ) ビニル ] ベンゼン  
(有機エレクトロルミネッセンス材料)

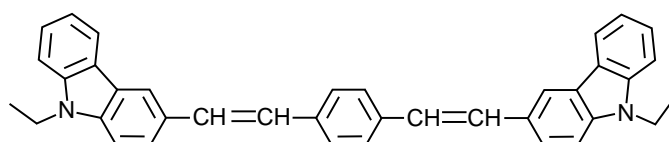

1,4-bis [2- (9-ethylcarbazol-3-yl) vinyl]  
benzene

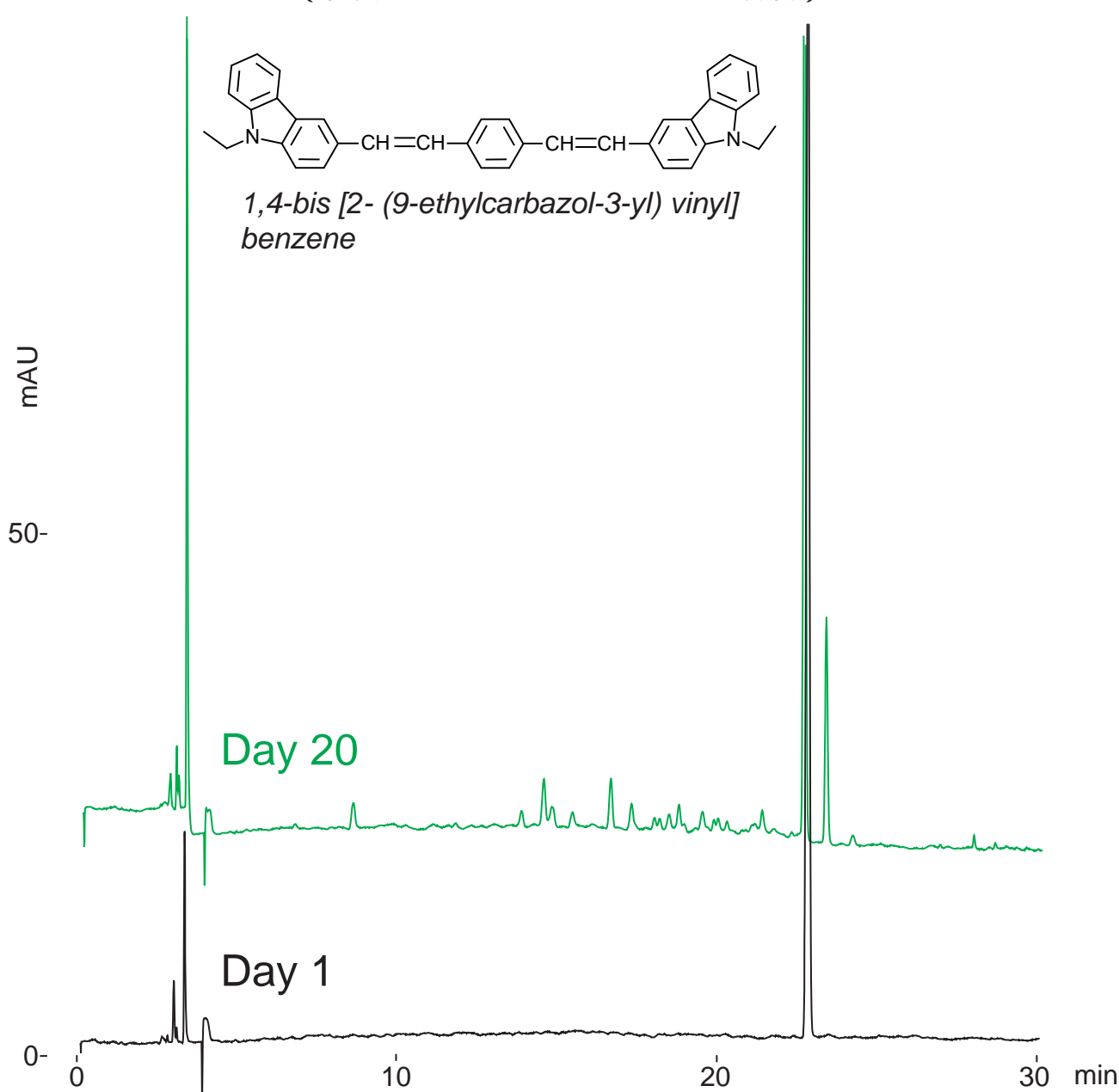

Unison UK-C8, 250 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-30min)

1.0mL/min (8MPa), 60°C, 260nm

10µL (0.05mg/mL, THF / formic acid = 100 / 1)